

| | |
|--------------|---------------|
| Pupil's name | Date of Birth |
|--------------|---------------|

Continuous Glucose Monitoring with 640G Insulin Pump

A continuous glucose monitor (CGM) is a device that measures interstitial glucose levels every 5 minutes and sends these readings to an insulin pump. As it is measuring interstitial glucose it can lag behind blood glucose levels. CGM does not take away the need for blood glucose testing but can be used to provide extra information about blood glucose trends.

The CGM is set with limits of acceptable glucose levels and will alert via the insulin pump if these limits are reached. These will be documented on the care plan clearly for staff reference.

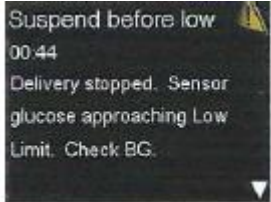

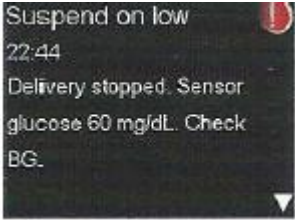
Continuous Glucose Monitoring (CGM)

| <u>Alert Message</u> | <u>What It Means</u> | <u>Action Required</u> |
|----------------------|--|--|
| Alert before low | <p>If Suspend before low is on, you will be alerted when the insulin is suspended.</p> <p>If Suspend before low is off, you will be alerted when the sensor predicts you will reach your low limit in 30 minutes</p> | <p>Do not treat your glucose based on the Sensor Glucose.</p> <p>Confirm using your blood glucose meter.</p> <p>Treat if necessary based on the hypo flowchart advice and continue to monitor.</p> |
| Alert on Low | Sensor glucose values are equal to or lower than your low limit. | |
| Alert before high | Sensor glucose reading is expected to reach the high glucose limit in the length of time that you set for the Time before high (5-30 minutes) | <p>Do not treat your glucose based on the Sensor Glucose.</p> <p>Confirm using your blood glucose meter.</p> |
| Alert on high | Sensor glucose values are equal to or higher than the high limit you set. | Treat if necessary based on the hyper flowchart advice and continue to monitor. |

- To acknowledge alerts/alarms follow on screen instructions.
- Do not calibrate sensor unless instructed to do so by parent.
- Blood glucose test required prior to bolus (snack, lunch and corrections).

| | |
|----------------------------|------|
| Signature of parent/carers | Date |
|----------------------------|------|

| | |
|--------------|---------------|
| Pupil's name | Date of Birth |
|--------------|---------------|

| SmartGuard suspend features | Impact on insulin delivery if suspend feature is turned on | Information displayed by the pump | |
|-----------------------------|---|---|--|
| Suspend before low | Insulin delivery is temporarily stopped if sensor glucose value is approaching your low limit | <p>You will receive this alert message and need to check your BG. Insulin delivery will remain suspended after the alert is cleared.</p>  <p>Of the alert is not cleared in 10 minutes, the pump will begin to siren.</p> | <p>After the alert or alarm message is cleared and insulin delivery has stopped, the home screen will display:</p> <ul style="list-style-type: none"> · Suspended before low or Suspended on low at the bottom of the screen · A shaded area to represent the time when the insulin has been suspended · A flashing SmartGuard technology icon.  |
| Suspend on low | Insulin delivery is temporarily stopped if sensor glucose value has reached or fallen below your low limit . | <p>You will receive this alarm message and need to check your BG. Insulin will remain suspended after the alarm is cleared.</p>  <p>If the alarm is not cleared after 10 minutes, the pump will begin to siren and an emergency message will appear on the pump screen.</p> | |

Signature of parent/carers

Date

| | |
|--------------|---------------|
| Pupil's name | Date of Birth |
|--------------|---------------|

Individual Patient Settings:

These are subject to change by the parent/carer and it is their responsibility to update the care plan when changes are made.

Low Limit: _____

High Limit: _____

Please tick box for what system is in use:

Smart Guard: ON ☐ OFF ☐

What the Sensor Icons Mean:

| | |
|------------|--|
| ↑ or ↓ | Sensor Glucose been rising or falling by 1 – 2 mmol/L over the last 20 minutes |
| ↑↑ or ↓↓ | Sensor Glucose been rising or falling by 2 – 3 mmol/L over the last 20 minutes |
| ↑↑↑ or ↓↓↓ | Sensor Glucose been rising or falling by 3 mmol/L over the last 20 minutes |

Optional SmartGuard Use:

If this function is switched on then the pump has the ability to:

- PREDICT when approaching low glucose levels 30 minutes in advance
- Automatically STOP insulin delivery before going hypo
- Automatically RESUME it when your glucose levels recover

Alert Message: If BG trending low or falls to low limit;

Suspend before low (Smart Guard use only)

Action:

- Allow pump to suspend and resume insulin delivery according to the BG levels.
- If bolus of insulin is required for food/correction and pump is suspended then the pump must be manually restarted before proceeding.
- Always do a manual BG check before resuming insulin and at the snack and mealtimes in insulin pump care plan.
- If hypo does occur then follow hypo flowchart as per the insulin pump care plan
- if the pump is suspended, resume 'basal deliver' once hypoglycaemia is resolved (blood glucose greater than 4.0mmol/L). To do this press 'suspend before low' and the pump will guide you through resuming the basal insulin.
- Calibrating sensor – always select 'no' when prompted.

| | |
|---------------------------|------|
| Signature of parent/carer | Date |
|---------------------------|------|